



FLOWASSIST 2-PIECE POUCH



HOW IT WORKS

Non-absorbent **FLOWASSIST** shields the skin and baseplate adhesive from exposure to stoma output.

Soft, flexible **FLOWASSIST** spout directs output away from the skin and into the pouch.

Baseplate can be rotated and positioned for the spout to align with stoma openings.

BENEFITS OF FLOWASSIST



- ✓ Leakage Protection
- ✓ Hydrocolloid Protection for Extended Wear Time
- ✓ Improved Skin Health
- ✓ Premium Security and Comfort

Designed to prevent problems before they begin.

Floating Flange

Allows easier pouch coupling with reduced pressure on the abdomen. Angled corners provide a secure grip for easy removal.

Radial Grooves

Comfortably conform to anatomical contours.

Contemporary Woven Fabric

Discreet, breathable, and water-repellent for everyday comfort.


Body-Contour Design

Follows natural curves for a secure fit with a low-rustle finish to minimize sound.



Belt Attachment Tabs on all Pouch types and Convex Baseplate for added security.

AVAILABLE IN

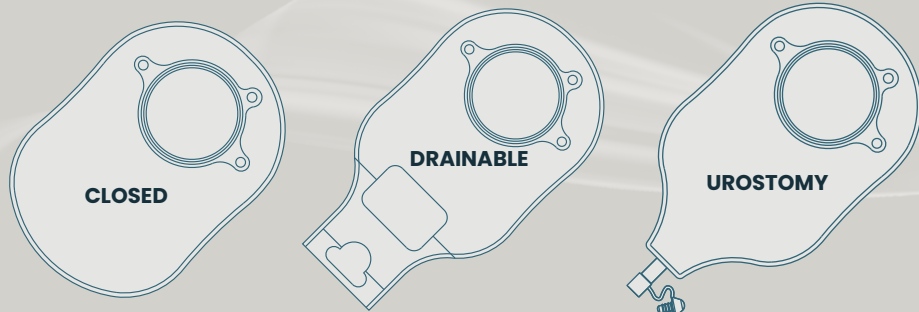

Closed and Drainable pouches are equipped with "Window" technology.

Closed pouches are fitted with an activated carbon filter that deodorizes gases while allowing them to escape unobstructed.

Drainable pouches include Moisture Protect Composition (MPC) filter technology - an additional protective layer that shields the carbon filter from liquid absorption.

Urostomy pouches feature a soft plug closure for easy emptying.

SIZING TABLE

HCPCS Code	Description	Size	Cut to Fit (Stoma Size)	SKU
A4410	FLOWASSIST Flat Baseplate	SMALL	19-25 mm / ¾ - 1"	OFA-FB57-S
		MEDIUM	25-32 mm / 1 - 1¼"	OFA-FB57-M
		LARGE	32-40 mm / 1¼ - 1 ⁹ / ₁₆ "	OFA-FB57-L
A4408	FLOWASSIST Convex Baseplate	SMALL	19-25 mm / ¾ - 1"	OFA-CB57-S
		MEDIUM	25-32mm / 1 - 1¼"	OFA-CB57-M
		LARGE	32-40mm / 1¼ - 1 ⁹ / ₁₆ "	OFA-CB57-L
A4425	FLOWASSIST Drainable 2-Piece Pouch	•	•	OFA-DP57
A4432	FLOWASSIST Urostomy 2-Piece Pouch	•	•	OFA-UP57
A4419	FLOWASSIST Closed 2-Piece Pouch	•	•	OFA-CP57
A4367	Ostoform Belt	120cm	•	OFA-OB-120



Watch Fitting Video

INTERESTED IN OUR PRODUCT?

Contact Ostoform's sales team
e: sales@ostofom.com



www.ostofom.com